



## Spore Methods Project – Plate Interpretation Guide


The purpose of this guide is to provide examples for how to enumerate colonies for dairy powder spore testing. In particular, we have included examples of spreading and swarming, which is common for spore-forming bacteria. This guide follows SMEDP<sup>1</sup> plate count interpretation for enumerating colonies on spread plate. As such, for plates with >25% spreading, estimate the count and mark “SPR” in the “Notes” column in the results spreadsheet.

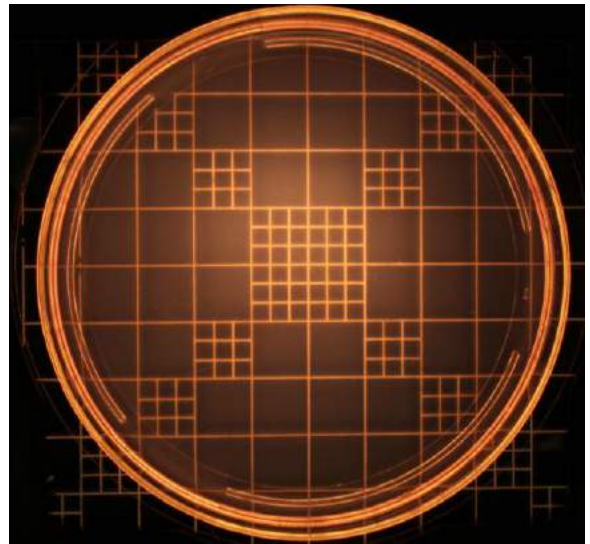

CFU = 0; no spreading

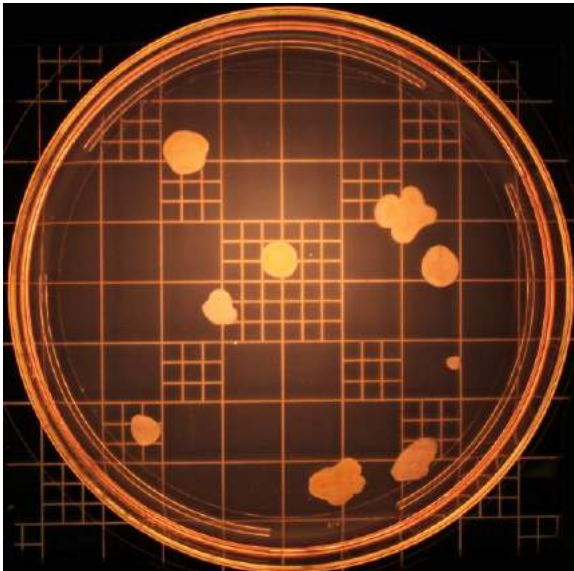

CFU = 9; no spreading

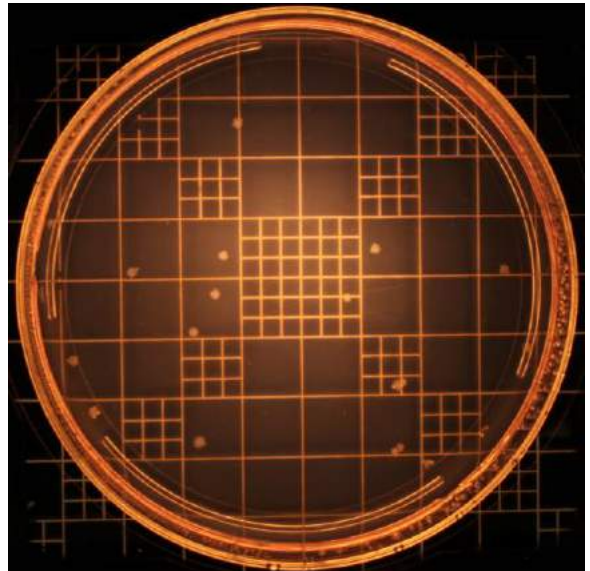

CFU = 16; no spreading




CFU = 3; ≤25% spreading




CFU = 1; ≤25% spreading




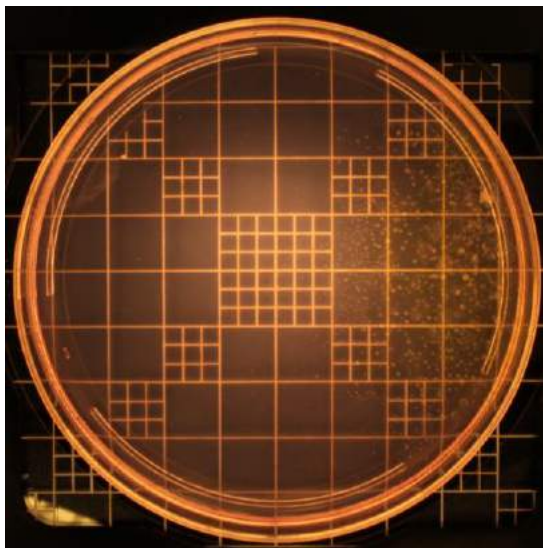


CFU = 1;  $\leq 25\%$  spreading

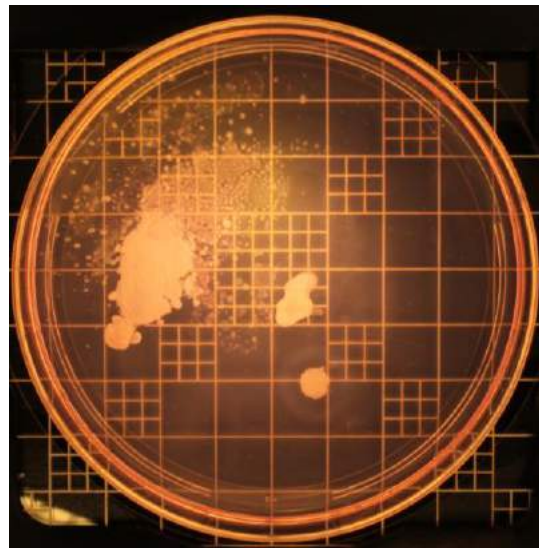


CFU = 3;  $\leq 25\%$  spreading

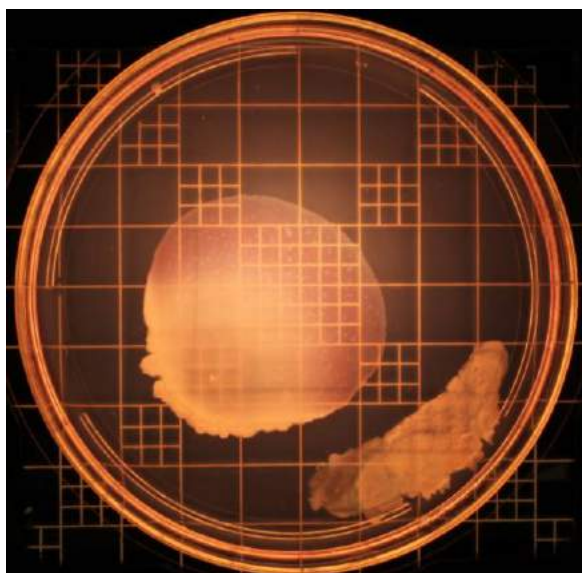


CFU = 2;  $> 25\%$  spreading = SPR

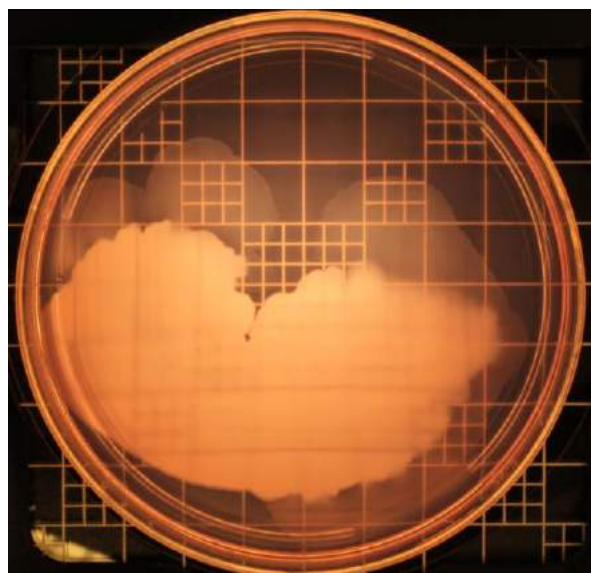

CFU = 1;  $> 25\%$  spreading = SPR

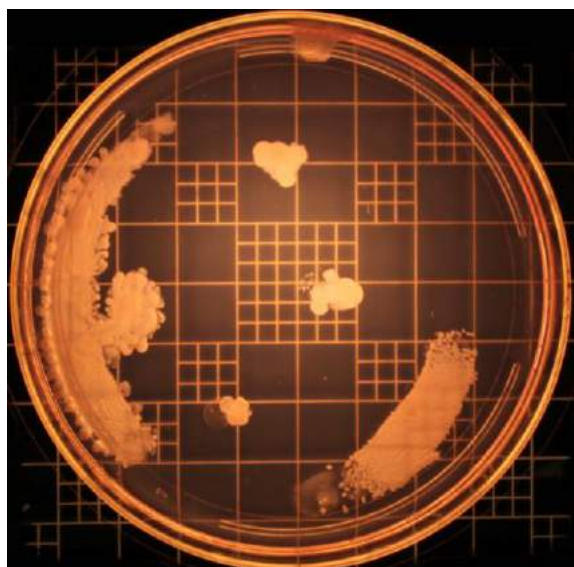

CFU = 6;  $\leq 25\%$  spreading

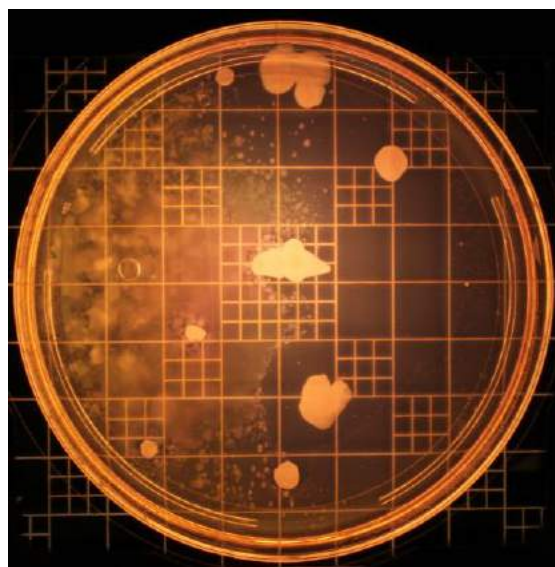

CFU = 9;  $> 25\%$  spreading = SPR

1. Laird, D. T., Gambrel-Lenarz, S. A., Scher, F. M., Graham, T. E. & Reddy, R. Microbiological count methods. in *Standard Methods for the Examination of Dairy Products* (eds. Wehr, H. & Frank, J.) 153–187 (American Public Health Association, 2004).